

Zephyr

Vascular Compression Device

Enhanced patient comfort while achieving patent haemostasis

The Zephyr vascular compression device has either a single or dual compression balloon that is easily adjustable to help facilitate patent haemostasis in the radial artery.

Simplify patent haemostasis with a clear radial balloon which provides excellent visibility for easy positioning and monitoring.



Partners in Healthcare Associated Blood Stream Infection Prevention Products and Strategies

Comfortable, convenient, cost effective

Zephyr Vascular Compression Devices


- The Zephyr 9100 device has a single compression balloon that is easily adjustable to help facilitate patent haemostasis in the radial artery. This technique has been shown to help avoid post-procedure chronic radial artery occlusion.¹
- The Zephyr 9200 device has two easily inflatable compression balloons, one for the radial artery and one for the ulnar artery. By compressing the ulnar artery along with radial artery, the Zephyr dual balloon reduces distal resistance to radial blood flow. This can help enhance antegrade radial blood flow during vascular compression to facilitate patent haemostasis, which has been shown to help avoid post-procedure chronic radial artery occlusion.¹

The Zephyr Vascular Compression Device may also be used in conjunction with other radial access products such as mBrace Wrist Support Device, WristAssist Wrist Support Device and Elevate Arm Board.

For use by medical professionals to promote haemostasis following a catheterisation or other puncture into a blood vessel in a patient's arm or leg, including: radial, brachial, dorsalis pedis, or tibial blood vessels; arterial or venous line or sheath removal; Haemodialysis, and in anticoagulation therapy.

Single or Dual Options

- The Zephyr helps clinicians achieve patent haemostasis with the clear balloon ensuring excellent puncture site visibility for easy positioning and monitoring

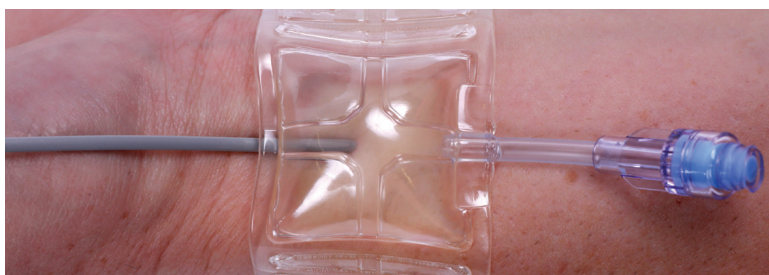

Easy One-handed Inflation and Deflation

- Every Zephyr Compression Device comes with a snap fit Syringe which allows for secure, controlled one handed inflation and deflation. The Zephyr valve is also compatible with a standard 20mL Luer Syringe



Enhanced patient comfort

- Simple and easy to use hook and loop closure on compliant elastomeric band for patient comfort giving firm compression for patent haemostasis available in multiple size options



References

1. Data on file

Product Codes

30/190-0100-00	Zephyr Single Balloon, Small, 21cm	50 / Carton
30/190-0100-10	Zephyr Single Balloon, Regular, 25cm	50 / Carton
30/190-0100-20	Zephyr Single Balloon, Large, 30cm	50 / Carton
30/190-0200-00	Zephyr Dual Balloon, Small, 21cm	50 / Carton
30/190-0200-10	Zephyr Dual Balloon, Regular, 25cm	50 / Carton
30/190-0200-20	Zephyr Dual Balloon, Large, 30cm	50 / Carton

Product Information



www.itlbiomedical.com/symbolglossary